

Features

A Perfect Gift for You & Your Loved Ones



Features



Heater

Press to turn on/off the heat function while massage continues.

Power

Press the POWER button to turn on the massager with heat. Auto timer will shut it off after 15 minutes. Press the Power button anytime to turn the unit off.


DR(Direction)

Press the DR button to change and maintain direction of the massage roller. When the unit is turned on, the massage roller will automatically change its rolling direction every 60 seconds.

Special features

- Special design Rolling massage system offers deep tissue Shiatsu massage.
- Flex handles let you easily target the right spot and apply pressure to the massage with your hands.
- Easy to operate. Equipped with Overheat Protection Device & programmed with auto timer to ensure safety.
- Come with a removable cotton cover for softer skin touch.
- Quiet operation, no noise. Durable leather-like cover.


Massage System

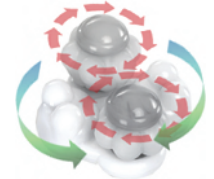

Silicone



Twisting

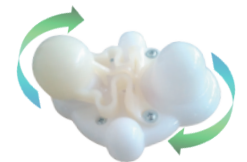

Springer

Full Body Massage

Full Body Massage

※ Please refer to the pictures below for a full body massage.



1 Neck and shoulder massage



2 Neck massage



3 Lower back massage



7 Back massage



8 Foot massage



9 Thigh massage



4 Abdomen massage



5 Hip massage



6 Calf massage



10 Ankle massage



11 Arm massage



12 Control the massage strength by adjusting the distance between armrests

Cautions and Notice

⚠ Cautions:

1. Do not use on an infant or invalid or unconscious person.
2. Do not use the product on body wounds or injuries.
3. Do not disassemble or modify the product on your own.
4. Do not drop the product or apply force on it as this may cause damage to the parts.
5. Do not use the products while sleeping.
6. For safety reason, do not use the product while driving.
7. Do not use the product continuously for more than 20 minutes.
8. The massager may cause uncomfortableness to patients with diabetes, acute illness, malignant tumor, high fever, skin disease, etc. Consult your doctor before using the product.
9. Keep close supervision when the product is used by, on or near Children.
10. Please do not wash the product or operate the product in moist environment, as it may cause electrical shock.
11. Please do not leave the product unattended when it is plugged in.
12. Please examine the covering before each use. And discard the product if any place shows any signs of deterioration.

⚠ Notice:

1. Confirm the voltage specified on the adaptor before using it.
2. Start the product after making sure the correct connection to transformer.
3. When using on neck or abdomen, always use it with wearing on clothes in order to protect your body skin.
4. Please remember to wipe the massager with a soft, damp cloth instead of abrasive cleaners.
5. Place the massager in a safe, cool and dry place when not in use.
6. Keep away from all solvents and harsh detergents.

Notice and Specification

Trouble Shooting:

1. When the massager is unable to turn the power on, that may be cause by poor insertion of power supply cord. Please ensure the power cord is properly plugged into the power socket.
2. When the massager functions are not working, that may be caused by overheating due to long period of using. Please be noted that the product will stop working automatically after continuous usage to prevent it from overheating. Please wait for more than one hour and the product could be used normally.

Warranty:

We warrant this product against defects in materials and /or workmanship under normal use for a period of ONE year from the date of purchase by the original purchaser

If the problem occurs under normal use and within the duration period of manufacturer' s limited warranty, the products is repaired or replaced with a new one by the manufacturer at its own options free of charge to the customers.

Specification :

Item Name: Neck and Shoulder Massager

Rated Voltage: DC 12V 2A

Power Consumption: 24W

Cover Material: Poly Urethane

Motion: Shiatsu with Heating

Packaging dimension: 40*16*27.5 CM